


CellSeal®

Closed-System Cryogenic Vial

Clear Access

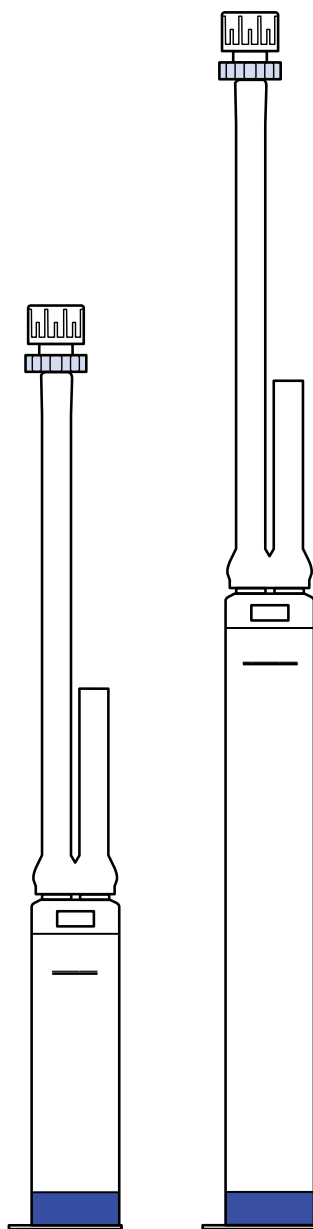


CELLSEAL® CLEAR ACCESS CLOSED-SYSTEM CRYOGENIC VIAL

NOTE: Handle the device using aseptic technique.

• RECOMMENDED MATERIALS

- Standard Luer connector syringe (10 mL syringe for 5 mL vial, and 5 mL syringe for 2 mL vial recommended)
- Exam gloves
- Processing rack
- Compatible sealing device (SEXTON RF Sealing System recommended)
- Clean scissors
- Alcohol swabs
- Compatible thawing accessory such as water bath or thawing system (SEXTON Thawing System recommended)
- Compatible cryogenic storage accessory, such as a box, cane, or controlled rate freezing container (SEXTON Cryogenic Box, Rack, and Freezing Container recommended)
- Cryogenic gloves

1. Wearing exam gloves, open the peelable pouches per user's validated methods. Remove the vial from the packaging. Place vial upright with tubing and port at the top in processing rack. Do not remove the vials from the package by the fill port. Handling product by the fill port weakens the tubing and fill port joint.

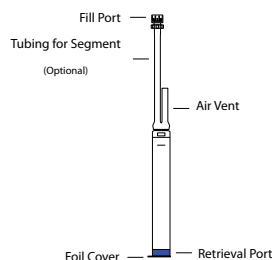

2. Working inside a biosafety hood or cleanroom, aseptically prepare the fill port by uncapping the Luer fitting. **NOTE:** Avoid pressing down on the cap when removing, as this can cause the cap to permanently stick to the fill port. It is also recommended to hold the stem of the Luer fitting while removing its cap to avoid weakening the tubing and fill port joint.


3. For loading the sample, use an appropriately sized syringe and needle to draw up the sample to be stored.

- a. If using a 5 mL vial, a 10 mL syringe is recommended.
- b. If using a 2 mL vial, a 5 mL syringe is recommended.


It is recommended to draw a small volume of air into the syringe first and then draw the sample. The air will be used to level the column of liquid in the tubing or to assist in pushing the entire sample into the vial. **Do not fill the vial with more than the stated maximum volume.**

4. Connect the syringe containing your sample to the fill port by turning the syringe clockwise. Once connected, push the sample into the vial. Do not use any non-standard Luer connectors; they may not fit with the fill port and cause spillage of the sample. To avoid dislodging the air vent filter, be sure to fill the vial at a moderate speed using the recommended needle and syringe size.

5. If a segment is desired, pull back on the syringe plunger to create a clear space in the tubing column above and below the liquid for the segment. Push the air into the tubing to place the liquid column at the desired level, or push the entire sample into the vial.

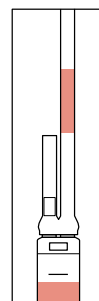

6. Using a sealer designed for use with ethylene-vinyl acetate tubing, seal the tubing as described in 6a or 6b below. **NOTE:** It is recommended to use SEXTON RF Sealing System for sealing of CellSeal vial's tubing.
 - a. If a segment is desired, place a seal above and below the liquid in the tubing. Place an additional seal below the bottom seal for optimal folding over of segment for storage.
 - b. If no segment is desired, place a seal near the vial body alone.
7. Seal the air vent tubing above the filter, making this a closed system. Once a closed system is created, use clean scissors to cut off any additional tubing. **NOTE:** The vial is intended to be frozen in the upright position only, and stored at cryogenic temperatures (as low as -196 degrees Celsius)

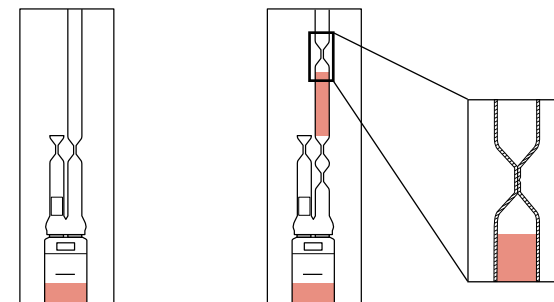

Seal Configuration for No Segment Seal Configuration for Segment Seal Close-Up View

8. Segment can be folded over for easy storage. Place vial and folded integral segment into box, cane, or controlled rate freezing container for storage. **NOTE:** SEXTON Cryogenic Storage Boxes, Racks, and Freezing Container are recommended.
9. When segment is needed, remove vial from cold storage and immediately cut off segment with clean scissors and thaw per user's validated sample thawing instructions.
10. When sample is needed, remove vial from cold storage and thaw per user's validated sample thawing instructions. **NOTE:** SEXTON Thawing System is recommended. **NOTE:** Vial is for single use only.
11. Cut open the vent tube with clean scissors.
12. Remove the foil on the bottom of the retrieval port. Thoroughly swab the retrieval port septum with a sterile alcohol wipe. Allow to air dry up to 1 minute before accessing.
13. Use a non-coring needle and/or proper anti-coring technique (as described below) to puncture the retrieval port septum and extract the sample. It is recommended to use an 18G non-coring needle for this port.

Proper Anti-Coring Technique:

 - A. Hold the needle with the opening of the needle tip away from the retrieval port septum.
 - B. Insert the needle into the septum at a 45-60° angle.
 - C. Increase the angle of the needle gradually as the needle enters the vial.
14. If used with infectious or hazardous bodily fluids, dispose of vial according to institutional guidelines for medical waste.

SEXTON

MANUFACTURER
Sexton Biotechnologies

1102 Indiana Avenue
Indianapolis, IN 46202
MADE IN USA



X = Quantity per packaging configuration

www.sextonbio.com

Revision: 2

Page 1 of 1